



*Evolution of Innovation!*

[www.thebetteroxygenmask.com](http://www.thebetteroxygenmask.com)



## One Oxy<sub>2</sub>Mask™ replaces multiple traditional oxygen therapy devices

Clinician efficiencies

Inventory reduction

Medical waste reduction

Oxygen consumption reduction

### Patient Safety

- Reduces the risk of CO<sub>2</sub> rebreathing
- Reduces the risk of emesis aspiration
- For nose or mouth breather

### Patient Comfort

- No build up of heat or humidity
- Reduces feelings of claustrophobia
- Unrestricted communication
- Access for fluids and oral medication



# Oxy<sub>2</sub>Mask™



## Changing the game with velocity vortices

- The updated Oxy<sub>2</sub>Mask design attaches the mushroom assembly directly to the diffuser wall as a baffle, with strategically positioned aerodynamic ports that act to direct the organized oxygen vortices with increased velocity specifically to the patient's nose and mouth

## Reliable performance


- The acclaimed flow dynamics of Oxy<sub>2</sub>Mask have evolved, continuing to offer the full range of FiO<sub>2</sub> delivery from 24% to 90%

## Dependable North American Supply Chain


- The updated design includes a snap-fit assembly and locally sourced materials, resulting in consistent product availability

## Specifications

- Delivers Oxygen flows from 1 to 15+ litres per minute - flush, 24% to 90% FiO<sub>2</sub> (Tyke size 1/4L to 5+, 22% to 65% FiO<sub>2</sub>)
- Is not made with natural rubber latex or phthalates such as DEHP
- 360° swivel adapter
- 25 masks/case
- Four sizes available:
  - Adult
  - Plus (25% larger mask)
  - Kid (3 - 10 years old, 33 - 70 lbs / 15 - 32 kg)
  - Tyke 6 months - 3 years old, 16-33 lbs / 7 - 15 kg)

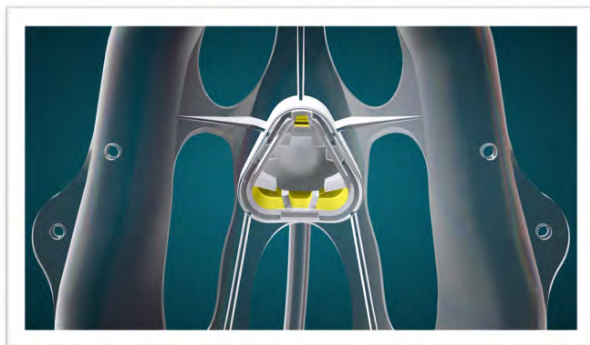

Oxy<sub>2</sub>Mask with enclosed diffuser and directional ports

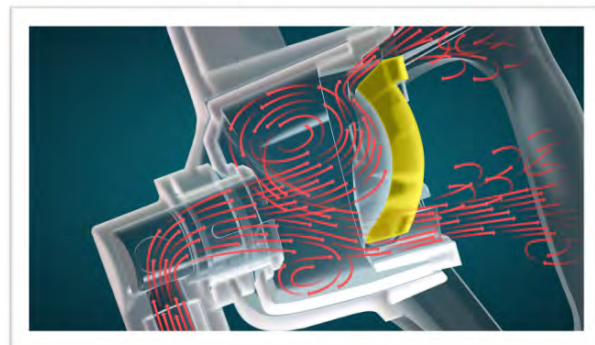

Side profile of the oxygen flow dynamics through the directional ports towards the mouth and nose (baffle highlighted).

Code	Description	Mask Size	Tubing Length	Elastic
OM-1125-8	Oxy <sub>2</sub> Mask™ Adult 7' tubing	Adult	7'	Single
OP-1125-8	Oxy <sub>2</sub> Mask™ Plus 7' tubing	Plus	7'	Single
OK-1125-8	Oxy <sub>2</sub> Mask™ Kid 7' tubing	Kid	7'	Single
OT-1025-8	Oxy <sub>2</sub> Mask™ Tyke 7' tubing	Tyke	7'	Single

*\*other variations, and non-PVC options available*



Respiratory  
Surgical  
Blades & Scalpels  
Dynamic Tissue Systems

**Manufactured by:**  
Southmedic Inc.  
50 Alliance Blvd.,  
Barrie, ON Canada L4M 5K3  
[www.thebetteroxygenmask.com](http://www.thebetteroxygenmask.com)

email: [custserv@southmedic.com](mailto:custserv@southmedic.com)  
1-800-463-7146 Toll Free in North America  
1-705-726-9383

ISO 13485

PO261 Sept. 2023

