Evolution of Innovation!



www.thebetteroxygenmask.com



One Oxy₂Mask[™] replaces multiple traditional oxygen therapy devices

SKU standardization

Inventory reduction

Medical waste reduction

Oxygen consumption reduction

Patient Safety

- Reduces the risk of CO2 rebreathing
- Reduces the risk of emesis aspiration
- For nose or mouth breather

Patient Comfort

- No build up of heat or humidity
- Reduces feelings of claustrophobia
- Unrestricted communication
- Access for fluids and oral medication



Oxy₂Mask[™]

Changing the game with velocity vortices

• The updated Oxy₂Mask design attaches the mushroom assembly directly to the diffuser wall as a baffle, with strategically positioned aerodynamic ports that act to direct the organized oxygen vortices with increased velocity specifically to the patient's nose and mouth

Reliable performance

- The acclaimed flow dynamics of Oxy_2Mask have evolved, continuing to offer the full range of FiO₂ delivery from 24% to 90%

Dependable North American Supply Chain


• The updated design includes a snap-fit assembly and locally sourced materials, resulting in consistent product availability

Specifications

- Delivers Oxygen flows from 1 to 15+ litres per minute flush, 24% to 90% FiO2 (Tyke size 1/4L to 5+, 22% to 65% FiO2)
- Is not made with natural rubber latex or phthalates such as DEHP
- 360° swivel adapter
- 25 masks/case
- Four sizes available:
 - Adult
 - Plus (25% larger mask)
 - Kid (3 10 years old, 33 70 lbs / 15 32 kg)
 - Tyke 6 months 3 years old, 16-33 lbs / 7 15 kg)



Oxy₂Mask with enclosed diffuser and directional ports

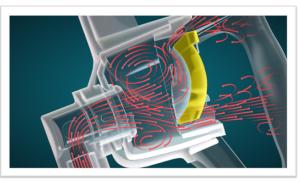

Side profile of the oxygen flow dynamics through the directional ports towards the mouth and nose (baffle highlighted).

Code	Description	Mask Size	Tubing Length	Elastic
OM-1125-8	Oxy₂Mask™ Adult 7' tubing	Adult	7′	Single
OP-1125-8	Oxy₂Mask [™] Plus 7' tubing	Plus	7'	Single
OK-1125-8	Oxy₂Mask [™] Kid 7' tubing	Kid	7'	Single
OT-1025-8	Oxy₂Mask [™] Tyke 7' tubing	Tyke	7'	Single

*other variations, and non-PVC options available





Manufactured by: Southmedic Inc. 50 Alliance Blvd., Barrie, ON Canada L4M 5K3 www.thebetteroxygenmask.com

email: custserv@southmedic.com 1-800-463-7146 Toll Free in North America 1-705-726-9383

ISO 13485

